

Adhesive Primer for Precious Metal

V-PRIMER

***Bond to Precious Metal
No more tin plating, heat treatment
For intraoral use***

What is V-PRIMER?

V-PRIMER is a one-component primer for bonding to precious metals (including Au/Ag/Pd alloys). It is based on the newly developed adhesive monomer called "VTD"*-a derivative of triazine-dithiol. Simple coat on precious alloy with V-PRIMER prior to application of Super-Bond C&B improves the durability of the bond. V-PRIMER eliminates the need for other bond-enhancing steps, such as heat treatment or tin plating.

* 6-(4-vinyl benzyl-n-propyl)amino-1,3,5 -triazine-2,4-dithiol

Made in Japan by




SUN MEDICAL CO., LTD.

571-2, Furutaka-cho, Moriyama, Shiga, 524-0044, Japan
Phone: 81-77-582-9981 Fax: 81-77-582-9984


Applications of V-PRIMER

V-PRIMER enhances the long-term bond strength of Super-Bond C&B to any gold or other precious metal restorations, such as:

1. Inlays, onlays, crowns, bridges
2. Adhesion bridges and adhesion splints
3. Precious metal posts
4. Orthodontic appliances to full metal crowns
5. Repair of worn acrylic veneers or fractured ceramometal crowns.


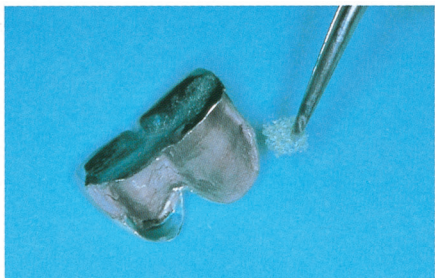


Application of V-PRIMER to crown

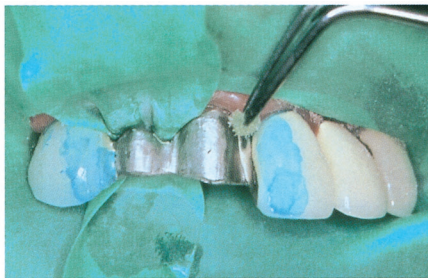

Intraoral application of V-PRIMER

Long Term Bond Strength (Shear strength)

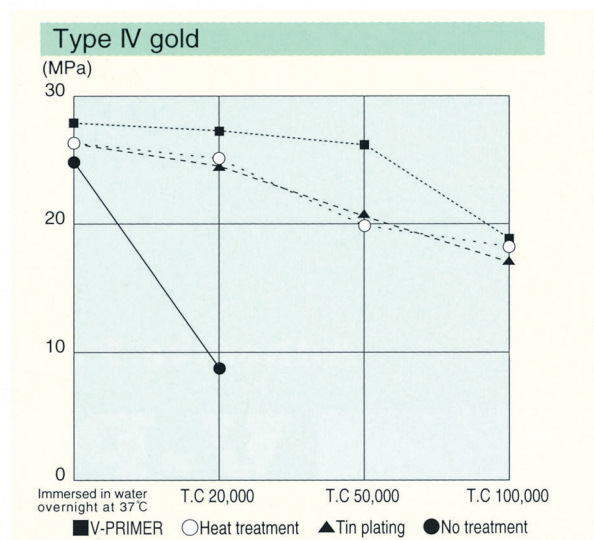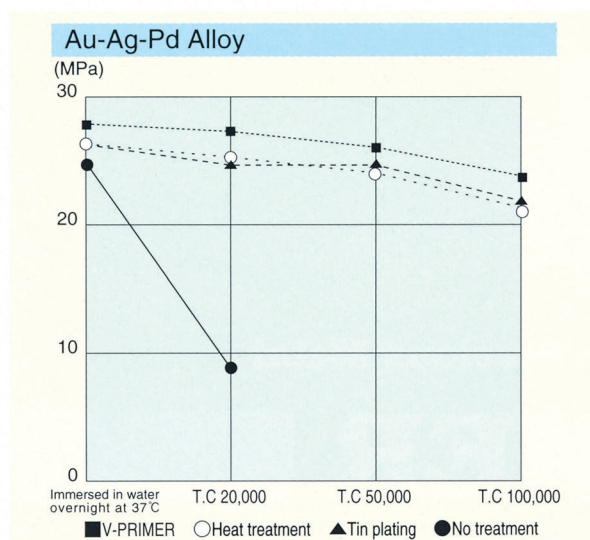

Contents



- Contents
 - V-PRIMER ----- 3mL
- Attached Documents
 - Instructions

Related product



Dental adhesive resin cement
Super-Bond
C&B

 **SUN MEDICAL CO.,LTD.**