## ■Related Products of Super-Bond C&B
(sold separately)

- **Quick Monomer**
  Optional monomer
- **Polymer (Clear, Esthetic, Opaque Ivory, Opaque Pink, L-Type Clear, L-Type Esthetic, L-Type Radiopaque)**
  Optional Polymer powders
- **Super-Bond Universal Ceramic Primer**
  Adhesive primer for ceramic/zirconia
- **V-PRIMER**
  Adhesive primer for precious metal alloys

## ■Accessories for Super-Bond C&B
(sold separately)

Dispensing Dish
Sponge (L･S)
Measuring Spoon (Standard) --------- for the standard mix ratio
Measuring Spoon (Small) ----- for 75% of the standard mix ratio
Measuring Spoon (Large) ---- for 120% of the standard mix ratio
Brush Handle (Straight)
Brush Handle (Bent)
Brush Tips (Blue) ---------- for Bulk-mix
Brush Tips (White･L) ---------- for Brush-dip
Brush Tips (White･S) ---------- for Brush-dip

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## CAUTION

①**Avoid contact**
Avoid contact with soft tissue, skin or eyes. A rubber dam is recommended for intraoral use. Dentist should use rubber or PVC dental gloves. Contaminated skin or mucosa should be wiped off immediately with alcohol and then thoroughly rinsed with running water, otherwise symptoms such as swelling may show. If Super-Bond C&B enters the eye, immediately rinse thoroughly with running water. The patient should be examined by an ophthalmologist. When the cured adhesive contacts soft tissue, polish the adhesive surface.
②**Be careful of acidity**
As the Red Activator and Green Activator are acid, avoid contact with soft tissue, skin or eyes, and do not allow patients to swallow them while applying and washing them with water.
③**Be careful of flammability**
Catalyst V and Monomer are flammable. Do not store where they may be exposed to open flame.
④**Clean spilled Catalyst V immediately with wet towels**
The Catalyst V reacts with oxygen. If absorbed by a flammable material, it may raise the temperature enough to cause smoldering. If the Catalyst V is spilled, wipe it up immediately with a **WET(not dry)** disposable towel. Then rinse the towel to kill the catalyst thoroughly in running water.

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## 1.What is Super-Bond C&B?

Super-Bond C&B is a self-cure dental adhesive resin cement based on MMA. It contains a high performance bonding monomer, "4-META", and a catalyst, "TBB". It shows excellent bond strength to tooth (enamel and dentin), metal*1, ceramic/zirconia*2 and dental resins.
For three decades it has been used extensively in a wide range of dental applications. During this long period it has earned an excellent reputation for pulpal safety.
Super-Bond C&B is widely known for its formation of a sound hybrid layer (resin impregnated layer) in both enamel and dentin. This layer reinforces the tooth surface against recurrent caries and prevents the post-operative hypersensitivity.

*1 Use V-PRIMER concurrently for precious metal alloys.
*2 Use Super-Bond Universal Ceramic Primer concurrently for ceramic/zirconia.

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## 2.Contents of the kit

| | |
|---|---|
| ①Monomer ---------- 10mL | ⑥Dispensing Dish ---------- 1 |
| ②Catalyst V ---------- 0.3mL | ⑦Measuring Spoon (Standard) ---------- 1 |
| ③Polymer (L-Type Clear)---------- 3g | ⑧Brush Handle (Bent)---------- 1 |
| ④Red Activator---------- 5mL | ⑨Brush Tips (White • L)---------- 10 |
| ⑤Green Activator---------- 5mL | ⑩Brush Tips (White • S) ---------- 10 |

■**Attached documents:**
Instructions / Pictorial Instruction Card

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# 3.Precautions

**Read all instructions thoroughly before use.**

## 3-1 Safety

**Please keep the following precautions for safe use.**
(Regarding Catalyst V, read 3-2 additionally.)

①**Applications**
Use Super-Bond C&B only for the applications recommended in this publication.

②**Past history of sensitivity**
Super-Bond C&B should not be used by clinicians or on patients who are methacrylic monomer-sensitive.

③**Symptomatic irritation**
Cease using Super-Bond C&B immediately, if signs of irritation such as rashes appear, and see a physician.

④**Avoid contact**
Avoid contact with soft tissue, skin or eyes. A rubber dam is recommended for intraoral use. Dentist should use rubber or PVC dental gloves. Contaminated skin or mucosa should be wiped off immediately with alcohol and then thoroughly rinsed with running water, otherwise symptoms such as swelling may appear. If Super-Bond C&B enters the eye, immediately rinse thoroughly with running water. The patient should be examined by an ophthalmologist. When the cured adhesive contacts soft tissue, polish the adhesive surface.

⑤**Be careful of acidity**
As the Red Activator and Green Activator are acid, avoid contact with soft tissue, skin or eyes, and do not allow patients swallow them during application or rinsing.

⑥**Pulp protection**
If the preparation approaches the pulp, apply a protective base.

⑦**Give care to flammability**
Catalyst V and Monomer are flammable. Do not store where they may be exposed to open flame.

## 3-2 Precautions on Catalyst V

**Catalyst V reacts with air and water to generate heat and lose activity. Please abide by the following.**

①**Storage conditions**
Avoid high temperature, high humidity and direct sunlight. The Catalyst should **NOT** be refrigerated. (The repeated temperature changes may shorten the Catalyst's shelf-life by causing the syringe to aspirate air.)
*After a long storage, the first drop of the Catalyst may be inactive, though the rest of the material remains active.
*The syringe is made of glass, therefore it must be handled with care to prevent shock, dropping, and other physical damage.

②**Cap closure**
The cap simply slides on and off. Recap the syringe immediately after each use. Air (oxygen and humidity) deactivates the Catalyst. Do not leave the cap off during the bonding procedure.

③**Screwing**
If the Catalyst does not come out of the syringe because of the tight screw, do not try too hard to trun it. The content may splash as the syringe breaks.

④**After use**
Unscrew the male-screw two turns counter-clockwise after each use to relieve pressure on the Catalyst. (Pressure buildup can cause leakage of the Catalyst or a crack of the syringe.)

⑤**Clean spilled Catalyst immediately with wet towels**
The Catalyst reacts with oxygen. If absorbed by a flammable material, it may raise the temperature enough to cause smoldering. If the Catalyst is spilled, wipe it up immediately with a **WET(not dry)** disposable towel. Then rinse the towel to kill the Catalyst thoroughly in running water.

⑥**Cleaning of the tip of the syringe**
Wipe the tip of the syringe with a dry gauze after each use to prevent residue buildup. Then rinse the gauze with water to kill any remaining activity. Buildup of the residue may prevent the cap from seating properly.

## 3-3 Storage

**Please take the following precautions to maintain the quality.**

①**Storage conditions**
As in the case of the Catalyst V, store the Monomer, Polymer, Red Activator and Green Activator in a cool, dark location. High temperature, high humidity and direct sunlight will shorten their shelf-life.

②**Volatility**
Monomer is highly volatile. Recap the bottle immediately.

③**Contamination**
Do not mix the bottle caps.

④**Dispensing dish cleaning**
After each use, wipe the dispensing dish with a solvent such as acetone before storing. (Alternatively, soak the dish in water. Then wash and dry.)

⑤**Brush cleaning**
The Brush Tips in the kit are for one-time use. Dispose of them after use. If any other multi-use brush is used, clean the brush immediately after each use with a solvent such as acetone before the resin has cured completely. Dry with paper tissue. If necessary, reshape the bristles so they will dry to a fine point.

## 3-4 Disposal

Dipose of empty containers, etc. in accordance with local regulations.

## 3-5 To Get the Best Results with Super-Bond C&B

①**Create and maintain a clean surface**
Oil, blood, saliva and biofilm will lower the bond strength. Clean the tooth and prosthesis thoroughly before cementing. After cleaning, take care to avoid re-contamination.

②**Dry the surfaces and prevent moisture contamination**
After cleaning, dry the surface adequately. A rubber dam is highly recommended, as it will reduce the chance of contamination by saliva, humid breath or blood.

③**Avoid eugenol-containing bases and cements**
Eugenol is a polymerization inhibitor. Therefore, eugenol-contained bases and cements should not be used with resin cements. To avoid cross-contamination, reserve a mixing dish exclusively for Super-Bond C&B. Do not use the same mixing dish for other adhesives.

④**Time constraints**
Super-Bond C&B's working and setting times are very different from those of traditional cements. Follow the instructions carefully to get the best results.

⑤**Do not re-use Polymer**
After using Super-Bond C&B in the Brush-dip technique, dispose of any excess Polymer left in the mixing dish. Do not return it into the container, as it has become contaminated with the Monomer.

⑥**Prosthesis Design**
To avoid stress concentration which encourages debonding regardless of the actual bond strength, design a prosthesis, such as wings of a bonded bridge, without thin unsupported area, which may flex during mastication.

***As in any dental treatment, the patient's individual constitution and the unique requirements of clinical case at hand must be considered before selecting materials and conditions for use.***

## 4-1 How to handle the containers, etc.

### Catalyst Syringe

**1. Decapping**
Pull off the Cap ① to open.


**2. Holding vertically**
Hold the Syringe vertical while expressing the catalyst to assure uniform drop size.

**3. Screwing**
Turn the screw ⑤ clockwise to express the catalyst.

**4. Relieving pressure**
Turn the screw ⑤ counter-clockwise two turns.

**5. Recapping**
Wipe the tip of the syringe with a dry gauze and replace the Cap ① immediately.
*See 3-2 Precautions on Catalyst V

### The structure of the Catalyst Syringe

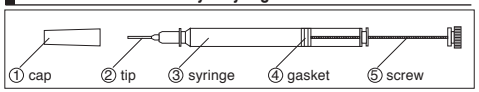■Plastic capsule over the refill Catalyst Syringe
The refill Catalyst Syringe comes in a plastic capsule to prevent breakage during transportation. Before use, please remove and dispose of the capsule.

### Monomer Bottle

**1. Decapping**
Unscrew the cap.

**2. Holding vertically**
Hold the bottle vertical while expressing the drops.

**3. Squeezing**
Squeeze the bottle to express drops.

**4. Recapping**
Recap immediately after use.
*See 3.3 ② Precaution against volatility.

### Polymer Jar

**1. Decapping**
Unscrew the cap.

**2. Opening aluminum seal**
Open the aluminum inner seal. (Only at initial use)

**3. Scooping powder**
Using the supplied spoon, scoop the Polymer powder.

**4. Leveling scoop**
Level the Polymer in the scoop using the cross-bar inside the jar.

**5. Recapping**
Replace the cap.

## 4-2 Operation Steps

**Super-Bond C&B can be used either with the Bulk-mix technique or the Brush-dip technique. Choose the appropriate technique by referring to the table below.**

### Comparison of Techniques

| | Bulk-mix technique | Brush-dip technique |
|---|---|---|
| **Outline of technique** | Polymer powder is mixed directly to the activated liquid*. | The powder/liquid ball is formed at the tip of a brush by dipping the tip first into the activated liquid* and then touching the Polymer powder. |
| **Comparison of the two techniques** | Use the powder/liquid mixture immediately. | The activated liquid must be used up within 5 minutes. |
| | Applicable to comparatively wide area. | Applicable to comparatively narrow area only. |
| | As the powder/liquid ratio is lower than that of Brush-dip technique, the working time is comparatively long but the curing is slow. | As the powder/liquid ratio is higher than that of Bulk-mix technique, the working time of mixed ball is comparatively short and the curing is fast. |

*Mixture of 4 drops of Monomer and 1 drop of Catalyst V

## Surface Preparation

It is essential that all surfaces to be bonded with Super-Bond C&B should be properly prepared. Preparation varies depending on the nature of the materials.

### Tooth Surface

**Moisture control**
Isolation by rubber dam or cotton roll is recommended.
*Remove scales beforehand

**Cleaning**
Remove contaminants and stains using a polishing brush and oil-free, fluoride-free pumice. Rinse thoroughly and dry.

**Surface treatment**
Apply the appropriate Activator using a sponge pledget or brush.

**Treatment time**

| Activator | Green | Red |
|---|---|---|
| Dentin | 5-10 sec | – |
| Enamel | – | 30 sec |

Do not use Red Activator on dentin.
Alternatively, enamel may be prepared with Green Activator for 30-60 seconds.

**Washing and drying**
Rinse thoroughly with water and dry.

### Ceramic/Zirconia Surface

**Cleaning**
Remove contaminants and stains from the ceramic/zirconia surface. Depending on the case at hand, cleaning may require abrasive wheels and points, or a prophy cup and fluoride-free, oil-free pumice.
(Applying the Red Activator will also help remove surface contaminants.)

**Washing and drying**
Rinse thoroughly and dry.
(Use of an ultrasonic cleaner is also effective.)

**Application of Super-Bond Universal Ceramic Primer**
Following the instructions that come with the primer, apply a single coat of Super-Bond Universal Ceramic Primer.

### Metal Surface

**Cleaning**
Remove contaminants and stains. Depending on the case at hand, cleaning may require abrasive wheels and points, or a prophy cup and fluoride-free, oil-free pumice.

**Sandblasting**
Air-abrade the surface with 50$\mu$m aluminum oxide.
*Rinse thoroughly and dry.
*After blasting, take care to avoid re-contamination. If the surface has become contaminated by oil or grease, the metal may be cleaned in an ultrasonic cleaner.
*If sandblasting is precluded (as, for example, in an intraoral repair), roughen the surface with an abrasive point.

**Protection of adjacent surfaces**
To facilitate removing the excess cement, apply a metal separating agent to surfaces that should not be bonded ... or protect the surfaces by wrapping with film.

**Application of V-PRIMER**
When bonding to precious metal, brush a thin coat of V-PRIMER onto the surface before applying Super-Bond C&B.
V-PRIMER does not enhance the bond to base metal alloys.

## 4-2 Operation Steps (Contd.)

## BULK-MIX TECHNIQUE

### Cooling the Dispensing Dish

Keep the Dispensing Dish in the refrigerator so it will be chilled when you use it. The recommended temperature range of the dish is 10 - 16˚C.

*If condensation forms on the dish when you remove it from the refrigerator, dry the well using an air syringe.

### Surface preparation

Surface preparation procedure varies depending on the materials being bonded (tooth,metal,or porcelain.)
(See "Surface Preparation" in the previous page.)

### Preparation of the Activated Liquid

**Dispensing the Monomer**

Hold the Monomer bottle vertical, and dispense the appropriate number of drops into the chilled Dispensing Dish.

**Dispensing the Catalyst V**


Hold the Catalyst syringe vertical, and turn the screw to dispense the proper number of drops to the Monomer. Stir lightly with a brush. This mixture is called "Activated Liquid".

### Mixing the Polymer

Using the supplied measuring spoon, add the Polymer powder to the Activated Liquid. Stir lightly with a brush.

### Application of the adhesive

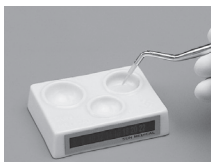Immediately after mixing, use a brush to apply the cement to the surface being bonded.
To deliver and apply Super-Bond more easily, recommend using Brush Tip (Blue) sold separately.

### Seating the restoration

Insert the restoration immediately. After confirming that it is completely seated, hold in position until the cement sets.

*Curing time is 8~10 min. at 37˚C for standard Polymer/Monomer ratio. The time varies with temperature, the Polymer type and the Polymer/Monomer ratio.

### Post treatment

Remove the excess cement. To facilitate this, protect beforehand unbonded surface properly and remove the excess resin timely.

Clean the Dispensing Dish used.

### Mixing ratio

| Monomer | Catalyst V | Polymer |
|---|---|---|
| 4 drops | 1 drop | 1 small cup of Measuring Spoon |
| 8 drops | 2 drops | 1 large cup of Measuring Spoon |

The Measuring Spoon (Standard) gives the standard Polymer/Monomer ratio.
The working time and the curing time can be controlled by changing the Polymer/Monomer ratio.
To change the Polymer/Monomer ratio, use an appropriate Measuring Spoon, which sizes are shown in the Table on the right.

### Relative volumes

1 corresponds to 0.2mL

| Measuring Spoon | Small cup | Large cup |
|---|---|---|
| Small | 0.75 | 1.5 |
| Standard | 1 | 2 |
| Large | 1.2 | 2.4 |

To increase the working time of a normal type of Polymer, use the Measuring Spoon (Small) sold separately.
To shorten the curing time using L-Type Polymer, use a Measuring Spoon (Large). The Measuring Spoon (Large) is sold separately.

### Key points to achieve good seating

**Work quickly, and seat the restoration before the mixture begins to gel.**

① Keep the Dispensing Dish chilled in the refrigerator. The recommended temperature: 10-16˚C.

② Mix Super-Bond C&B at the last moment before bonding.

③ In case you want to increase the working time further, use an L-Type Polymer and/or reduce the Polymer/Monomer ratio to 0.75 using a Measuring Spoon (Small) sold separately.

④ In case you want to shorten the curing time while assuring the reasonable working time, use an L-Type Polymer and increase the Polymer/Monomer ratio to 1.2.
A Measuring Spoon (Large) with large cups (1.2 / 2.4) is available separately.

## BRUSH-DIP TECHNIQUE

### Dispensing the Polymer

Dispense an appropriate amount of Polymer powder into the well "P" of the Dispensing Dish.

### Surface preparation

Surface preparation procedures vary depending on the materials being bonded (tooth, metal, or porcelain.)
(See 4-2 "Surface Preparation".)


### Preparation of the Activated Liquid

**Dispensing the Monomer**

Hold the Monomer bottle vertical and dispense the appropriate number of drops into well "L" of the Dispensing Dish.

**Dispensing the Catalyst V**

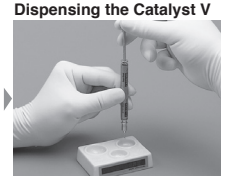Hold the Catalyst syringe vertical, and turn the screw to dispense the proper number of drops into the Monomer. Stir lightly with a brush. This mixture is called "Activated Liquid".

### Mixing ratio

| Monomer | Catalyst V |
|---|---|
| 4 drops | 1 drop |
| 8 drops | 2 drops |


### Application of the Activated Liquid

Brush the liquid onto the surface to be bonded.

*The Activated Liquid decomposes gradually and loses activity. Use it within 5 min. after preparation.

### Brush-dip Procedure

**Dipping the brush**

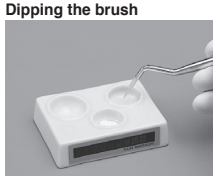Dip the Brush Tip into the Activated Liquid in the well "L". Eliminate excessive liquid by touching the brush tip to the edge of the well.

*When you repeat the procedure, clean up the brush with gauze before you dip it.

**Forming the ball**

Touch the brush to the Polymer powder in the well "P". A small ball of powder will be picked up on the wet tip of the brush.

**Applying the ball**

Brush the powder ball onto the prewet surface being bonded. As soon as it touches the surface, the powder will spread out to create a creamy, homogeneous layer. If necessary, repeat the procedure until the entire surface is covered with the cement.

### Seating the restoration

Insert the restoration immediately. After confirming that it is completely seated, hold in position until the cement sets.

*Curing time is 5~6 min. at 37˚C.

### Post treatment

Remove the excess cement.
To facilitate this, protect any adjacent surfaces you will not be bonding, and remove the excess resin before it sets.
Clean the Dispensing Dish used.

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